

Canis Major™

WET DENTAL /
TREATMENT
LIFT TABLE

Matrx
by MIDMARK

Matrx by Midmark is pleased to introduce the newest Canis Major model for the treatment room. The same adjustable-height lift table is now available with wet treatment tub. The table is designed to eliminate lifting to the treatment table and to allow for a variety of operator positions. The Canis Major wet table adjusts to individual needs whether your staff is short, tall, works sitting or standing.

The tabletop lowers to within 7" of the floor and raises to 44" high. The shallow wet tub is excellent for all treatment and dental procedures.

The conveniently located temperature controlled ultralight 8' sprayer extends beyond the table length for complete coverage and easy cleaning. The removable hair trap and easy clean drain make maintenance a breeze. Connections provided for easy installation to hot and cold water and drain.

Specifications:

The Canis Major is manufactured to exacting standards of fit and finish.

The shallow tub table and rack is constructed of stainless steel. Other metal components are powder coated for both chip and scratch resistance.

The back modesty cover fully encases all moving parts.

Temperature controlled water with ultralight 8' sprayer.

Removable hair trap and easy clean drain.

Two (2) GFI electric outlets.

Internal sensor and peristaltic pump for automatic draining.

300 lb Maximum Lift

Powered by medical-grade electric actuator for maintenance-free quiet, smooth lift

Lowers to within 7" of floor and rises to 44"

Hand operated motion controls

110 AC @ 10 AMP power

Overall dimensions:

45" Model: 62"L x 60"H x 25"W

60" Model: 77"L x 60"H x 25"W

**NEW LIFTING WET TABLE
EXCELLENT FOR DENTALS**
NOW AVAILABLE WITH SCALE!



Available in 60"
or 45" Tabletops

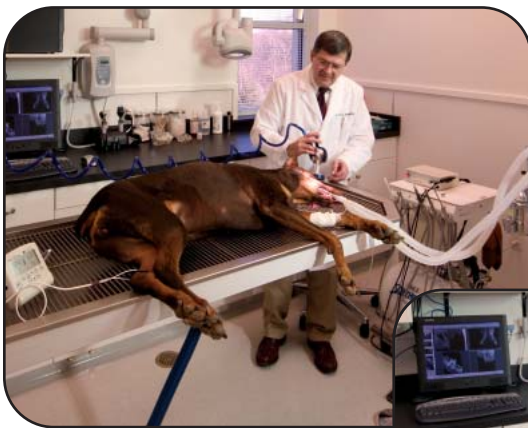


Photo: Cummings Veterinary Hospital,
Clifford Cummings, DVM
Easton, PA



Side View

Matrx
by MIDMARK

1-800-MIDMARK | midmark.com